

interscience



BagMixer® 400

The world's best selling lab blender

Why choose **BagMixer®** 400!

Best-seller

- Sold in more than 80 countries
- 30 years experience

High-quality

- Optimal analysis with the sterile **BagFilter®**
- In compliance with ISO 7218 standard
- Made in France by **interscience**

Brushless motor

- No wearing parts
- - 50% power consumption
- Extra silent

Safety

- The handle closes the blender and starts the blending
- **Q-Tight®**: no leaks from the top of the bag
- Security drip tray

PATENTED

Click & Clean® removable paddles

- Remove the paddles and easily clean the blending chamber



Easy cleaning

- All stainless steel
- Door opens at 270°
- Removable **Click & Clean®** paddles

Guarantee

- Shock absorbers and window door **life-time** guarantee
- 3 year guarantee (after registration)

PATENTED


Adjustable blending power

- Up to 28 kg of pressure: the hardest samples don't resist! e.g.: nuts, pills, seeds, dog biscuits etc.



Optimal extraction

- From as little as 30 seconds
- Easy control with the digital screen (speed, time, power)

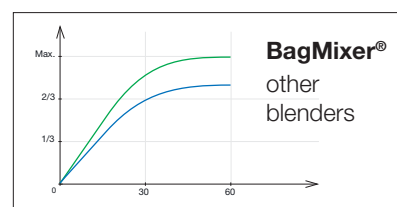

More than 10,000 **BagMixer®** in use every day!

Made in France by **interscience**,
the only manufacturer of blenders and sterile bags for microbiological analysis

Some features are limited to specific models.
For complete technical info, please check online:
www.interscience.fr



4 BagMixer® 400 models:

The **BagMixer®** lab blender guarantees optimal extraction by the action of its two paddles. Adapted to any sample (food, agricultural, chemical, medical, environmental) to be analyzed between 50 ml and 400 ml, it guarantees homogenization in 30 to 60 seconds.



BagMixer® 400 P

> The essential

Ref. 021 230



BagMixer® 400 W

> Window door

Ref. 022 230



BEST SELLER

BagMixer® 400 VW

> Adjustable speed

Ref. 023 230



BagMixer® 400 CC

> Click & Clean®

Ref. 024 230

Technical specifications at the back of the leaflet →

BagMixer® packs:

Starter Pack



1 BagMixer®

+ Sterile bags: 500 BagFilter®
+ 1,000 BagLight®

+ Accessories: 50 BagClip®
+ 1 BagOpen® + 1 BagRack®

Classic Pack



1 BagMixer®

+ Sterile bags: 500 BagFilter®
+ 1,000 BagLight® + 500 BagPage®

+ Accessories: 50 BagClip®
+ 1 BagOpen® + 1 BagRack®
+ 1 PetriPile® + 1 BagPipet®
+ 1,000 BagTips® Regular
+ 6 Clinet® 101

Premium Pack



1 BagMixer® & 1 BabyGravimat®

+ Sterile bags: 1,000 BagFilter®
+ 2,000 BagLight® + 1,000 BagPage®

+ Accessories: 50 BagClip®
+ 1 BagOpen® + 2 BagRack®
+ 2 PetriPile® + 1 BagPipet®
+ 1,000 BagTips® Regular
+ 6 Clinet® 101

BagSystem® complete range for sample preparation before analysis



BagFilter®
Collect



Gravimat®
Dilute



BagMixer®
Blend



BagTools®
Pipet

Technical specifications

		BagMixer® 400 P	BagMixer® 400 W	BagMixer® 400 VW	BagMixer® 400 CC®
Reference		021 230	022 230	023 230	024 230
Key features	Removable and autoclavable paddles, Click & Clean® system	-	-	-	✓
	Blending power indicator	-	-	-	✓
	Adjustable blending power (up to 28 kg of pressure)	-	-	✓	✓
	Digital screen	-	-	✓	✓
	Leak free Q-tight® system	-	-	✓	✓
	Side-by-side paddle stop	-	-	✓	✓
	Silent brushless motor	-	-	✓	✓
	Window door	-	✓	✓	✓
	Security drip tray	in option	✓	✓	✓
	270° opening door	✓	✓	✓	✓
	All stainless steel	✓	✓	✓	✓
	Compact and ergonomic	✓	✓	✓	✓
Technical specifications	Door weight	1050 g	950 g		
	Blending time	30, 60, 90, 120, 150, 180, 210 sec., ∞		from 1 sec. to 59 min, ∞	
	Blending speed	fixed speed: 8 strokes/sec.		variable speed: 4, 6, 8, 10 strokes/sec.	
	Dimension (L x l x h)	39 x 26 x 29 cm		42 x 26 x 26.5 cm	
	Weight	16.5 kg			
	Power	~100-127V / ~220-240V - 50-60Hz		~100-240V - 50-60Hz	
	Capacity	50-400 ml			
Guarantees	Made in France	✓	✓	✓	✓
	Life-time guarantee	Unbreakable shock absorbers and window door			
	3 year guarantee (after registration)	✓	✓	✓	✓
Manufactured in compliance with					 

Discover our complete range for microbiology



BagFilter®
Sterile bags



Gravimat®
Gravimetric dilutors



BagMixer®
Lab blenders



easy**Spiral**[®]
Spiral platers



Scan®
Colony counters

Your local distributor:

interscience INTERNATIONAL

30, Ch. Bois des Arpents - 78860 St Nom - FRANCE
T: +33 (0)1 34 62 62 61 - F: +33 (0)1 34 62 43 03
info@interscience.fr - www.interscience.fr

interscience USA & CANADA

199 Weymouth ST - ROCKLAND - MA 02370 - USA
T: 781 792 2133 - F: 781 792 2134
info@intersciencelab.com - www.intersciencelab.com